WEAVE Open-Time Phase 2. Instructions for Applicants

Tiny link: https://a.ing.iac.es/s/?21dfdf

- · Important links and documentation
- Communication with the PI
- Summary of actions
- LIFU mode: catalogue preparation
 - Target declination and right ascension limits
 - Total requested time
 - Preparation of the catalogue
 - Filling out the catalogue template with targets
 - TARGPROG
 - TARGID and TARGNAME
 - PROGTEMP
 - Instrument configuration
 - OB length
 - Arm exposures
 - Binning in the spectral direction
 - Duplicate this OB 'x' times
 - OBSTEMP
 - DIMM seeing
 - Transparency
 - Elevation angle
 - Moon distance
 - Sky brightness
 - GAIA specific columns
 - IFU_PA_REQUEST
 - IFU_DITHER
 - Summary of table columns and examples
 - Aladin overlays of the LIFU footprints
 - Catalogue upload to WASP
- LIFU mode: OB preparation
 - Download the validated catalogue
 - Generation of OBs
 - OB upload to WASP
- Observations and follow up

WEAVE Open-Time Phase 2. Instructions for Applicants

Important links and documentation

WEAVE Automated Submission Platform (WASP): http://wasp.ast.cam.ac.uk/

IFU workflow package: https://ingbitbucket.ing.iac.es/projects/WVSWG/repos/ifu/browse/workflow

WEAVE Open Time programme: https://www.ing.iac.es/astronomy/observing/WEAVE-OT.html

A detailed description of how the observations are processed and eventually released can be found in the WEAVE data model documentation:

- 1. Login to WASP or the OR.
- 2. Click on ("Tools" and select "Data Model").
- In section "Resources", click on the "Documentation" tab and download WEAVE-ICD-030: "Interface between SPA and SWG/QAG (info for WASP etc)".

Communication with the PI

The WEAVE Automated Submission Platform (WASP) automatically emails the principal investigators (PIs) with submission deadlines, validation acknowledgements and any other important information concerning the preparation of phase 2. The first two emails contain the following (the subject line is provided):

1. An email containing the credentials to log onto WASP:

```
Subject: [WEAVE] Account created for WEAVE observation preparation
```

2. A following email contains information about the opening of the OB submission window:

Subject: [WEAVE-WASP] 20XXXX WASP upload window now open

General assistance can be provided at the following email:

weave_open_time_support@ing.iac.es

Summary of actions

GREEN: Actions by the PI

PURPLE: Actions by the PI if familiarised with the package "IFU Workflow", or by ING on behalf of the PI, if assistance is requested by the PI to weave_o pen_time_support@ing.iac.es

Step	Action
1	Download of catalogue template from WASP (optional)
2	Fill out the catalogue template with targets following instructions on this page
3	Process the populated catalogue using the IFU workflow
4	Upload the processed catalogue to WASP
5	Download the validated catalogue from WASP
6	Generate the XML files (OBs) using the IFU workflow
7	Upload OBs to WASP
8	Select OBs on WASP for observing

The instructions below are mainly aimed at PIs who will request assistance to:

weave_open_time_support@ing.iac.es

LIFU mode: catalogue preparation

Target declination and right ascension limits

The current target declination limit is -25°. Objects with a declination lower than that cannot be observed.

WASP sets limits to allowed target RA and Dec, and it will not validate your OB if your target coordinates lie outside the limits for that trimester. See below the allowed RAs for each trimester.

Trimester A1	
2h (30°) 20h (300°)	
Expanded range: 2h-3h-4h18h-19h-20h 30°-45°-60°270°-285°-30)0°
Trimester A2	
9h (135°) 2h (30°)	
Expanded range: 9h-10h-11h0h-1h-2h 135°-150°-165°0°-15°-30)0
Trimester B1	
15h (225°) 9h (135°)	
Expanded RA range: 15h-16h-17h7h-8h-9h 225°-240°-255°105°-120	°-135°

Trimester B2

20h (300°) 16h (240°) Expanded RA range: 20h-21h-22h-...-14h-15h-16h 300°-315°-330°-...-210°-225°-240°


Total requested time

The total length (counted as science exposure time) of all the submitted OBs can't exceed the science time allocation of your proposal (see the time allocations).

Preparation of the catalogue

The catalogue can be prepared in three alternative ways:

1. Filling in the FITS catalogue template provided on WASP. Log on to WASP with your credentials. Go to "My survey" and click on "Download template".

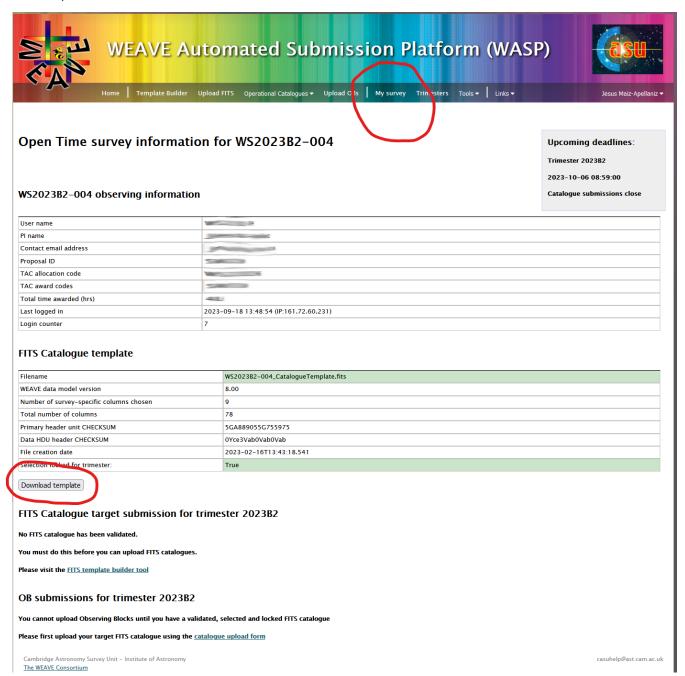

- 2. Converting the above FITS catalogue template into CSV and filling in.
- 3. Filling in a simplified CSV version of the catalogue template. Below it's an example of this simplified template and 4 catalogue rows that each describe an OB definition:

```
<filename: WS2023B2-001_target_list.csv>

TARGID,TARGNAME,GAIA_RA,GAIA_DEC,PROGTEMP,OBSTEMP,IFU_DITHER
IC768_1,IC768,76.7836980,18.1489034,61331.3,DAFEA,3
GHMon,GHMon_NW,32.591451,23.910124,60881,IADEE,4
SN2024gtp,SN2024gtp,102.644083,+68.252231,43331.3,DACEA,6
M101HIIregionC1,M101HIIregion,161.2533346,-11.7653604,51001.7+,
AAAEA,0
```

Filling out the catalogue template with targets

Any of the three catalogue templates above (FITS, CSV, simplified CSV) can be filled in with targets. These are the requested columns (the rest are optional). If not specified, the data release version and epoch of the GAIA coordinates will be DR3 and 2016.0 respectively.

TARGPROG (optional)
TARGID
TARGNAME
GAIA_ID (if available)
GAIA_DR (optional)
GAIA_RA
GAIA_DEC
GAIA_EPOCH (optional)
GAIA_PMRA (if available)
GAIA_PMDEC (if available)
GAIA_PARAL (if available)
PROGTEMP
OBSTEMP
IFU_DITHER
IFU_PA_REQUEST (if necessary)


TARGPROG

TARGPROG is an optional column, to be filled out at the discretion of the catalogue creator.

TARGID and TARGNAME


This parameter is used to group IFU observations of the same target, in cases where stacks are required. This helps the Core Processing System (CPS) responsible for obtaining the data, to identify cases where the same astrophysical target is observed but the OBs executed were not related (via for example the "chained" directive).

An example would be LIFU observations of the core of M33. If a user requires 3 OBs, each with different dither positions, then the CPS could not ordinarily stack these data, because they do not share the same Central CNAME (CNAME is the WEAVE object name from coordinates). In the input FITS catalogue, these observations are represented by 5,427 rows: 603 fibres, dithered 3 times for 3 OBs. Each of these rows must be tied together by a common TARGNAME, e.g. "M33 bulge". This indicates to the CPS that these observations should be evaluated for stacking.

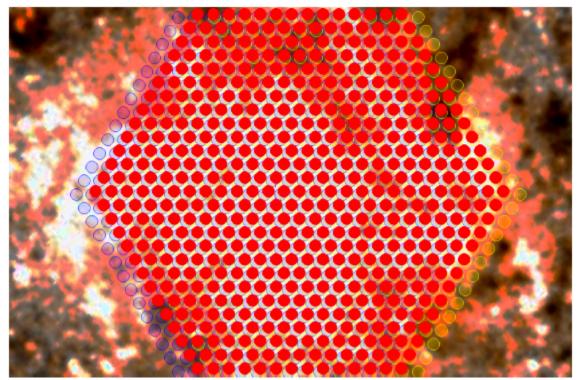

Above: Three LIFU dither positions centred on a putative target. This information must be encoded in the input FITS catalogue, but grouped by common TARGNAME.

It is not always true that IFU observations with common TARGNAME will be stacked. Only common TARGNAMEs with sufficient overlap will be stacked by the CPS. However, if Contributed Data Products exist to create larger mosaics from these data, then they should use the common TARGNAME to group L1 products.

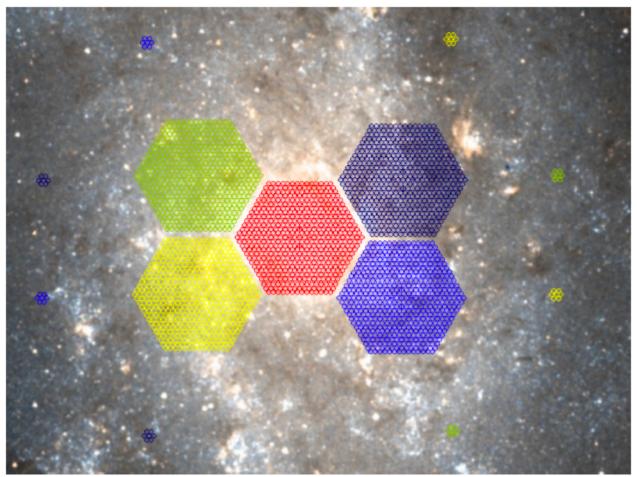

Above: Example LIFU coverage of an extended source. Each colour here represents an OB. Within each pointing there would be a series of dithers (as per previous Figure). Whilst these pointings might share the same TARGNAME ("M33"), the CPS would recognise these as a mosaic and not stack them into a single LIFU data cube.

The TARGID is the target identifier and it is mandatory that this column is filled in.

For IFU fields, this is the OB-specific descriptor of the field. In the first example above, this could be a simple numerical identifier for each OB:

```
OB1: TARGNAME = "M33bulge" TARGID = "M33bulge1"

OB2: TARGNAME = "M33bulge" TARGID = "M33bulge2"

OB3: TARGNAME = "M33bulge" TARGID = "M33bulge3"
```

For the second example, TARGID could be more descriptive:

```
OB1: TARGNAME = "M33" TARGID = "M33 bulge"
OB2: TARGNAME = "M33" TARGID = "M33 disc NE"
OB3: TARGNAME = "M33" TARGID = "M33 disc NW"
```

PROGTEMP

The PROGTEMP code is an integral part of describing how the instrument is configured. This parameter encodes the requested instrument configuration, OB length, exposure time, spectral binning, cloning requirements and probabilistic connection between these clones. We refer registered WASP users to the PROGTEMP form: http://wasp.ast.cam.ac.uk/progtemp (there is also a copy at ING: https://a.ing.iac.es/weave/progtemp.php).

Please use this exposure time calculator for your signal-to-noise calculations: https://a.ing.iac.es/signalWEAVE/. There is more information, including a downloadable version, at WEAVE Exposure Time Calculator.

Instrument configuration

Check first what instrument modes are available in the ING announcement of opportunity.

OB length

This PROGTEMP component specifies the overall OB length (inclusive of overheads). Please be aware that not all OB lengths are available.

Arm exposures

It defines how the time within the OB that you wish to observe the target is divided up, by specifying the number of exposures and exposure time within the stipulated OB length. The listed options of exposure time splits correspond to science open-shutter times (note that every exposure incur a 3-minute overhead).

IMPORTANT: It is not recommended unlocking the arms. Note that different exposure splits (unlocking the arms) in the red and the blue arms will have two consequences: firstly, no dithering will be possible, and secondly, OB length will increase as the observing system cannot execute different exposure splits simultaneously but one only exposure in both arms at a time. For example, requesting one exposure in red, and 12 in blue would cause a significant mismatch in overheads between the two arms.

Binning in the spectral direction

Only listed values are allowed, all other values are forbidden. As an aside, note that there is no facility to request spatial binning. Such a mode would compromise the spectral extraction performance and crosstalk characterisation within the L1 pipeline.

Duplicate this OB 'x' times

This component is entirely optional – validation checks will not fail if there is no 'x' value in the PROGTEMP. If an applicant wishes to duplicate the observations of this target, then duplication would result in the OB this target lives in being duplicated 'x' times, or in other words, one only OB will be generated, but it will be cloned 'x' times when ingested into the observing system.

Furthermore, if PIs want to *chain* these OBs together, such that the observation of one of these clones increases the chance that the remaining OBs in this group are observed, then a "+" will be added to the end of PROGTEMP. The "+" instructs the WEAVE scheduler to increase the (internal, programme-specific) weight of all OBs in this group. This provides a mechanism to ensure that once a series of observations are started, completion of the chain becomes a progressively higher priority above other OBs from this programme.

OBSTEMP

Whilst PROGTEMP deals with "how" a target is observed, OBSTEMP deals with "when" a target is observed – namely setting the observational constraints required to optimally extract scientific information from the observation.

These constraints represent "worst-case" observing parameters – targets will be observed under these conditions or better. Note, however, that the worst case represents a scenario when all worst-case observational constraints meet simultaneously. Pls should note also that the probability of observing targets with highly restrictive OBSTEMP codes may be quite low.

We refer registered WASP users to the OBSTEMP form: http://wasp.ast.cam.ac.uk/obstemp (there is also a copy at ING: https://a.ing.iac.es/weave/obstemp.php).

DIMM seeing

Seeing is provided by a DIMM monitor located nearby the WHT, and it is measured as the full-width half-maximum (FWHM) brightness profile of a star if located at airmass 1.00. Note that IFU observations are more tolerant of poorer seeing (diameter of the fibres is 2.6"). An instrument-level limit of seeing <3" is adopted due to degradation of guiding quality in poorer seeing.

Transparency

This component defines the minimum transparency (T) required for the observation. In the absence of cloud and Saharan dust, the V-band transparency (i. e. the fraction of top-of-the-atmosphere light reaching the WHT) is \sim 0.88, and likely scales as (0.88)^airmass (i.e. T \sim 0.83 at airmass 1.5, T \sim 0.77 at airmass 2). Extinction due to Saharan dust rarely exceeds a few tenths of a magnitude (i.e. reduces transparency by $<\sim$ a few tenths) and varies little with position on the sky. Thin cloud, reducing transparency by a few tenths, can be patchy, so only an all-sky average can be predicted (based on the guidestar signal from previous observations). A limit on transparency (T) > 0.4) is adopted to avoid various performance degradations such as guidance on attenuated guide stars, delivery of science data with significantly reduced signal-to-noise ratio and increased risk of precipitation on the telescope optics.

Elevation angle

Users may set the minimum elevation angle (and thus airmass) that their targets should be observed with. An instrument-wide limit of airmass < 3 will be set due to several effects degrading the quality of data at higher airmass: increasing differential refraction; higher extinction; worse light pollution. Please remember to check for the visibility of your target using http://c.ing.iac.es/staralt/ when placing an elevation limit.

Moon distance

Users may specify the minimum angular distance between the moon and the target. Observations under low moon distance may suffer from sky background light gradients across the field, as well as higher levels of reflected light causing spurious signals within the focal plane.

An instrument-wide lower limit of 30 degrees ensures sky brightness gradients and scattered light effects are kept to a minimum. This limit applies to any observations set with "0 degrees (no constraint)".

Sky brightness

The maximum V-band surface brightness of the sky (mag / sq. arcsec) required to observe the target.

The dark-of-moon sky brightness varies by a few tenths of a mag depending on ecliptic and galactic latitude and phase of the solar cycle. With the moon up, the sky can be up to ~4 mag brighter than dark-of-moon, depending on lunar phase, elevation and angular distance from the target. Science observations are carried out mainly during astronomical night, but those OBs tolerating moonlit skies can also use the latest bit of evening twilight, or the earliest bit of morning twilight – the only relevant criterion is sky surface brightness. No observations will be carried out when the sky brightness exceeds that at full moon.

GAIA specific columns

We strongly urge PIs to use the DR3 reference epoch for WEAVE targets.


Objects that are not detected by Gaia must have their coordinates converted onto the Gaia Reference Frame DR3. It is the responsibility of the PI to verify that their provided coordinates are both consistent with the ICRS.

Parallax and proper motions are particularly important for construction of the WEAVE CNAME, to ensure that successive observations (with varying separation in time) of the same target are correctly attributed to the same CNAME.

Information from Gaia is expected to be provided exactly as it is retrieved from the Gaia archive, so parameter definitions, data formats and units remain the same. Please consult the Gaia webpages to further understand their data model.

IFU_PA_REQUEST

IFU_PA_REQUEST allows the catalogue provider to specify a rotation angle (position angle) of the LIFU, if this is required. Any dithering requested for the field is applied to the rotated LIFU frame. If a specific PA value is not requested, then the following defaults are applied: IFU_PA=0° (for declinations lower than 28.7°) and IFU_PA=180° (for declinations higher than 28.7°). Note that the requested IFU_PA could be changed by the IFU workflow software if a suitable guiding star cannot be found at the requested PA. Below it's a plot which shows how the PA is changed to find a suitable guiding star.

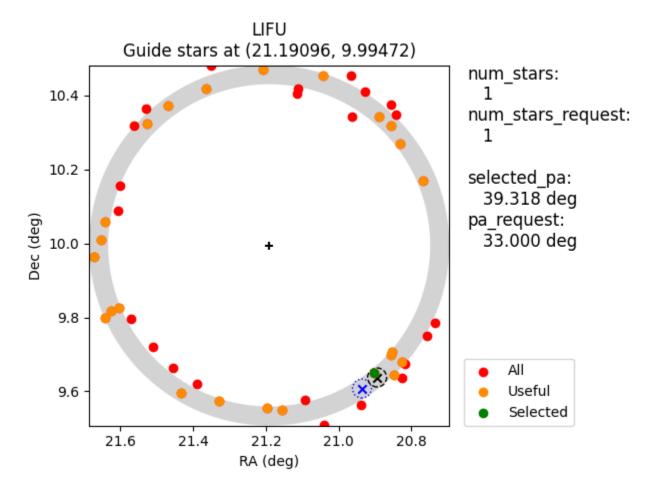

IMPORTANT: During the IFU observation preparation using the IFU workflow software, an analysis is performed on putative LIFU pointings to determine if a rotation is required and provides the optimal value IFU_PA should take. Rotation is generally used in cases where the default or requested PA results in no viable guide stars falling within the autoguider field of view, or when the rotation during an exposure is likely to exceed the angular limits of the rotator. This means that it's not always possible that the IFU workflow software can provide the IFU_PA_REQUEST value in the catalogue. The table below shows the allowed IFU_PA_REQUEST (=PA in the table) ranges PIs can safely apply for. Maximum value of PA is 360°.


	h_min=50.28 (X=1.3) OBSTEMP:A PA range (°)			5.58 (X=1.4) MP:B) h_min=41.81 (X=1.5) h_min=38.68 (X=1.6) OBSTEMP:C OBSTEMP:D					3.75 (X=1.8) MP:E	h_min=25.00 (X=2.4) OBSTEMP:F		
Dec (°)			PA range (°) PA range (°)		PA range (°)		PA range (°)		PA range (°)		PA range (°)		
-25									225 <pa<135< th=""><th>-135 0 +135</th><th>245<pa<115< th=""><th>-115 0 +115</th></pa<115<></th></pa<135<>	-135 0 +135	245 <pa<115< th=""><th>-115 0 +115</th></pa<115<>	-115 0 +115	
-10	220 <pa<140< td=""><td>-140 0 +140</td><td>235<pa<125< td=""><td>-125 0 +125</td><td>245<pa<115< td=""><td>-115 0 +115</td><td>245<pa<115< td=""><td>-115 0 +115</td><td>250<pa<110< td=""><td>-110 0 +110</td><td>260<pa<100< td=""><td>-100 0 +100</td></pa<100<></td></pa<110<></td></pa<115<></td></pa<115<></td></pa<125<></td></pa<140<>	-140 0 +140	235 <pa<125< td=""><td>-125 0 +125</td><td>245<pa<115< td=""><td>-115 0 +115</td><td>245<pa<115< td=""><td>-115 0 +115</td><td>250<pa<110< td=""><td>-110 0 +110</td><td>260<pa<100< td=""><td>-100 0 +100</td></pa<100<></td></pa<110<></td></pa<115<></td></pa<115<></td></pa<125<>	-125 0 +125	245 <pa<115< td=""><td>-115 0 +115</td><td>245<pa<115< td=""><td>-115 0 +115</td><td>250<pa<110< td=""><td>-110 0 +110</td><td>260<pa<100< td=""><td>-100 0 +100</td></pa<100<></td></pa<110<></td></pa<115<></td></pa<115<>	-115 0 +115	245 <pa<115< td=""><td>-115 0 +115</td><td>250<pa<110< td=""><td>-110 0 +110</td><td>260<pa<100< td=""><td>-100 0 +100</td></pa<100<></td></pa<110<></td></pa<115<>	-115 0 +115	250 <pa<110< td=""><td>-110 0 +110</td><td>260<pa<100< td=""><td>-100 0 +100</td></pa<100<></td></pa<110<>	-110 0 +110	260 <pa<100< td=""><td>-100 0 +100</td></pa<100<>	-100 0 +100	
0	250 <pa<110< td=""><td>-110 0 +110</td><td>255<pa<105< td=""><td>-105 0 +105</td><td>255<pa<105< td=""><td>-105 0 +105</td><td>260<pa<100< td=""><td>-100 0 +100</td><td>260<pa<100< td=""><td>-100 0 +100</td><td>270<pa<90< td=""><td>-90 0 +90</td></pa<90<></td></pa<100<></td></pa<100<></td></pa<105<></td></pa<105<></td></pa<110<>	-110 0 +110	255 <pa<105< td=""><td>-105 0 +105</td><td>255<pa<105< td=""><td>-105 0 +105</td><td>260<pa<100< td=""><td>-100 0 +100</td><td>260<pa<100< td=""><td>-100 0 +100</td><td>270<pa<90< td=""><td>-90 0 +90</td></pa<90<></td></pa<100<></td></pa<100<></td></pa<105<></td></pa<105<>	-105 0 +105	255 <pa<105< td=""><td>-105 0 +105</td><td>260<pa<100< td=""><td>-100 0 +100</td><td>260<pa<100< td=""><td>-100 0 +100</td><td>270<pa<90< td=""><td>-90 0 +90</td></pa<90<></td></pa<100<></td></pa<100<></td></pa<105<>	-105 0 +105	260 <pa<100< td=""><td>-100 0 +100</td><td>260<pa<100< td=""><td>-100 0 +100</td><td>270<pa<90< td=""><td>-90 0 +90</td></pa<90<></td></pa<100<></td></pa<100<>	-100 0 +100	260 <pa<100< td=""><td>-100 0 +100</td><td>270<pa<90< td=""><td>-90 0 +90</td></pa<90<></td></pa<100<>	-100 0 +100	270 <pa<90< td=""><td>-90 0 +90</td></pa<90<>	-90 0 +90	
+15	270 <pa<90< td=""><td>-90 0 +90</td><td>270<pa<90< td=""><td>-90 0 +90</td><td>270<pa<90< td=""><td>-90 0 +90</td><td>270<pa<90< td=""><td>-90 0 +90</td><td>275<pa<85< td=""><td>-85 0 +85</td><td>275<pa<85< td=""><td>-85 0 +85</td></pa<85<></td></pa<85<></td></pa<90<></td></pa<90<></td></pa<90<></td></pa<90<>	-90 0 +90	270 <pa<90< td=""><td>-90 0 +90</td><td>270<pa<90< td=""><td>-90 0 +90</td><td>270<pa<90< td=""><td>-90 0 +90</td><td>275<pa<85< td=""><td>-85 0 +85</td><td>275<pa<85< td=""><td>-85 0 +85</td></pa<85<></td></pa<85<></td></pa<90<></td></pa<90<></td></pa<90<>	-90 0 +90	270 <pa<90< td=""><td>-90 0 +90</td><td>270<pa<90< td=""><td>-90 0 +90</td><td>275<pa<85< td=""><td>-85 0 +85</td><td>275<pa<85< td=""><td>-85 0 +85</td></pa<85<></td></pa<85<></td></pa<90<></td></pa<90<>	-90 0 +90	270 <pa<90< td=""><td>-90 0 +90</td><td>275<pa<85< td=""><td>-85 0 +85</td><td>275<pa<85< td=""><td>-85 0 +85</td></pa<85<></td></pa<85<></td></pa<90<>	-90 0 +90	275 <pa<85< td=""><td>-85 0 +85</td><td>275<pa<85< td=""><td>-85 0 +85</td></pa<85<></td></pa<85<>	-85 0 +85	275 <pa<85< td=""><td>-85 0 +85</td></pa<85<>	-85 0 +85	
+28.70	295 <pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td></pa<65<></td></pa<65<></td></pa<65<></td></pa<65<></td></pa<65<></td></pa<65<>	-65 0 +65	295 <pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td></pa<65<></td></pa<65<></td></pa<65<></td></pa<65<></td></pa<65<>	-65 0 +65	295 <pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td></pa<65<></td></pa<65<></td></pa<65<></td></pa<65<>	-65 0 +65	295 <pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td></pa<65<></td></pa<65<></td></pa<65<>	-65 0 +65	295 <pa<65< td=""><td>-65 0 +65</td><td>295<pa<65< td=""><td>-65 0 +65</td></pa<65<></td></pa<65<>	-65 0 +65	295 <pa<65< td=""><td>-65 0 +65</td></pa<65<>	-65 0 +65	
+30	130 <pa<230< td=""><td>-130 180 +130</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<225<></td></pa<230<></td></pa<230<></td></pa<230<></td></pa<230<>	-130 180 +130	130 <pa<230< td=""><td>-130 180 +130</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<225<></td></pa<230<></td></pa<230<></td></pa<230<>	-130 180 +130	130 <pa<230< td=""><td>-130 180 +130</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<225<></td></pa<230<></td></pa<230<>	-130 180 +130	130 <pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<225<></td></pa<230<>	-130 180 +130	135 <pa<225< td=""><td>-135 180 +135</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<225<>	-135 180 +135	135 <pa<225< td=""><td>-135 180 +135</td></pa<225<>	-135 180 +135	
+40	115 <pa<245< td=""><td>-115 180 +115</td><td>120<pa<240< td=""><td>-120 180 +120</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<230<></td></pa<240<></td></pa<240<></td></pa<240<></td></pa<245<>	-115 180 +115	120 <pa<240< td=""><td>-120 180 +120</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<230<></td></pa<240<></td></pa<240<></td></pa<240<>	-120 180 +120	120 <pa<240< td=""><td>-125 180 +125</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<230<></td></pa<240<></td></pa<240<>	-125 180 +125	120 <pa<240< td=""><td>-125 180 +125</td><td>130<pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<230<></td></pa<240<>	-125 180 +125	130 <pa<230< td=""><td>-130 180 +130</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<230<>	-130 180 +130	135 <pa<225< td=""><td>-135 180 +135</td></pa<225<>	-135 180 +135	
+50	105 <pa<255< td=""><td>-105 180 +105</td><td>100<pa<260< td=""><td>-100 180 +100</td><td>115<pa<245< td=""><td>-115 180 +115</td><td>120<pa<240< td=""><td>-120 180 +120</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<240<></td></pa<240<></td></pa<245<></td></pa<260<></td></pa<255<>	-105 180 +105	100 <pa<260< td=""><td>-100 180 +100</td><td>115<pa<245< td=""><td>-115 180 +115</td><td>120<pa<240< td=""><td>-120 180 +120</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<240<></td></pa<240<></td></pa<245<></td></pa<260<>	-100 180 +100	115 <pa<245< td=""><td>-115 180 +115</td><td>120<pa<240< td=""><td>-120 180 +120</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<240<></td></pa<240<></td></pa<245<>	-115 180 +115	120 <pa<240< td=""><td>-120 180 +120</td><td>120<pa<240< td=""><td>-125 180 +125</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<240<></td></pa<240<>	-120 180 +120	120 <pa<240< td=""><td>-125 180 +125</td><td>135<pa<225< td=""><td>-135 180 +135</td></pa<225<></td></pa<240<>	-125 180 +125	135 <pa<225< td=""><td>-135 180 +135</td></pa<225<>	-135 180 +135	
+60	80 <pa<280< td=""><td>-80 180 +80</td><td>95<pa<265< td=""><td>-95 180 +95</td><td>105<pa<255< td=""><td>-105 180 +105</td><td>110<pa<250< td=""><td>-110 180 +110</td><td>120<pa<240< td=""><td>-120 180 +120</td><td>130<pa<230< td=""><td>-130 180 +130</td></pa<230<></td></pa<240<></td></pa<250<></td></pa<255<></td></pa<265<></td></pa<280<>	-80 180 +80	95 <pa<265< td=""><td>-95 180 +95</td><td>105<pa<255< td=""><td>-105 180 +105</td><td>110<pa<250< td=""><td>-110 180 +110</td><td>120<pa<240< td=""><td>-120 180 +120</td><td>130<pa<230< td=""><td>-130 180 +130</td></pa<230<></td></pa<240<></td></pa<250<></td></pa<255<></td></pa<265<>	-95 180 +95	105 <pa<255< td=""><td>-105 180 +105</td><td>110<pa<250< td=""><td>-110 180 +110</td><td>120<pa<240< td=""><td>-120 180 +120</td><td>130<pa<230< td=""><td>-130 180 +130</td></pa<230<></td></pa<240<></td></pa<250<></td></pa<255<>	-105 180 +105	110 <pa<250< td=""><td>-110 180 +110</td><td>120<pa<240< td=""><td>-120 180 +120</td><td>130<pa<230< td=""><td>-130 180 +130</td></pa<230<></td></pa<240<></td></pa<250<>	-110 180 +110	120 <pa<240< td=""><td>-120 180 +120</td><td>130<pa<230< td=""><td>-130 180 +130</td></pa<230<></td></pa<240<>	-120 180 +120	130 <pa<230< td=""><td>-130 180 +130</td></pa<230<>	-130 180 +130	
+70			65 <pa<295< td=""><td>-65 180 +65</td><td>85<pa<275< td=""><td>-85 180 +85</td><td>95<pa<265< td=""><td>-95 180 +95</td><td>110<pa<250< td=""><td>-110 180 +110</td><td>120<pa<240< td=""><td>-125 180 +125</td></pa<240<></td></pa<250<></td></pa<265<></td></pa<275<></td></pa<295<>	-65 180 +65	85 <pa<275< td=""><td>-85 180 +85</td><td>95<pa<265< td=""><td>-95 180 +95</td><td>110<pa<250< td=""><td>-110 180 +110</td><td>120<pa<240< td=""><td>-125 180 +125</td></pa<240<></td></pa<250<></td></pa<265<></td></pa<275<>	-85 180 +85	95 <pa<265< td=""><td>-95 180 +95</td><td>110<pa<250< td=""><td>-110 180 +110</td><td>120<pa<240< td=""><td>-125 180 +125</td></pa<240<></td></pa<250<></td></pa<265<>	-95 180 +95	110 <pa<250< td=""><td>-110 180 +110</td><td>120<pa<240< td=""><td>-125 180 +125</td></pa<240<></td></pa<250<>	-110 180 +110	120 <pa<240< td=""><td>-125 180 +125</td></pa<240<>	-125 180 +125	
+80							35 <pa<325< td=""><td>-35 180 +35</td><td>90<pa<270< td=""><td>-90 180 +90</td><td>140<pa<220< td=""><td>-140 180 +140</td></pa<220<></td></pa<270<></td></pa<325<>	-35 180 +35	90 <pa<270< td=""><td>-90 180 +90</td><td>140<pa<220< td=""><td>-140 180 +140</td></pa<220<></td></pa<270<>	-90 180 +90	140 <pa<220< td=""><td>-140 180 +140</td></pa<220<>	-140 180 +140	

IFU_DITHER

IFU_DITHER allows the PI to specify the dithering strategy for their observations, but may alternatively request not to dither their observations, with the understanding that this will not provide full spatial coverage of their requested field. PIs may request the preset -3, 3, 4, 5 and 6-dither patterns. The LIFU preset patterns will be properly rotated according to the position angle of the observation.

IMPORTANT: Exposure options for which fixed dither patterns are available in IFU modes are 3x, 4x, 5x, 6x.

Dither code	Description		owed odes Dither pattern (coordinates given in arcsec)										c)	
-1	Custom dither	LI	LIFU Defined by users via the RA and Dec of each position, as described in the OB <field> element (Section)</field>										n, as	
			LIFU	X					0.0	0000				
0	No dither	MOS	LIFU	у					0.0	0000				
	No didici	MOS	mIFU	X					0.0	0000				
			mire	у					0.0	0000				
				X	0	.000	0		-	1.9494	4		-0.	9747
3	3-dither preset	LI	FU	у	0	.000	0		(0.0000)		-1.	6882
		mIFU		X	0	0.0000			0.0000			-0.9208		
				у	0	.000	0	1.00			.0632		0.5316	
		LIFU		X	0.000	0	-1.94		4	-0.9747		-0.9747		0.9747
4	4-dither preset			У	0.000	0	0.00)	-1.6882		32 1.6882		1.6882
4	4-diffici preset	mIFU		X	0.0000		0.0000	0.0000 -0.920		.920			0.9208	
				у	0.000	0000 1.0		1.0632			6	. (0.5316	
		LIFU		X	0.0000)	-1.9494		-0.9	9747	-0.	.9747		0.9747
5	5-dither preset			У	0.0000)	0.0	000	-1.0	5882	1.	6882	$oxed{oxed}$	1.6882
,	5-diffici preset	mIFU		X	0.0000	$\overline{}$	0.0000					9208		0.9208
				У	0.0000			632	Ь——	316		5316	Ļ	-0.5316
		T.T	FU	X	0.0000		9494	-0.97	$\overline{}$	-0.97	\rightarrow	0.974	\rightarrow	-2.9240
6	6-dither (3 and			У	0.0000		0000	-1.68		1.68	$\overline{}$	1.688	$\overline{}$	1.6882
	-3 combined)	mI	FU	X	0.0000				208	0.92	$\overline{}$			0.0000
				У	0.0000		0632	0.53		0.53	16	-0.531	_	-1.0632
		LIFU		X	-0.9747		0.9747				-2.9240			
-3	3-dither alt	Life		У		1.6882		1.6882			\rightarrow	1.6882		
	preset	mI	mIFU			208		0.9208				0.0000		
					0.5	316			-0.5	316		-1.0632		

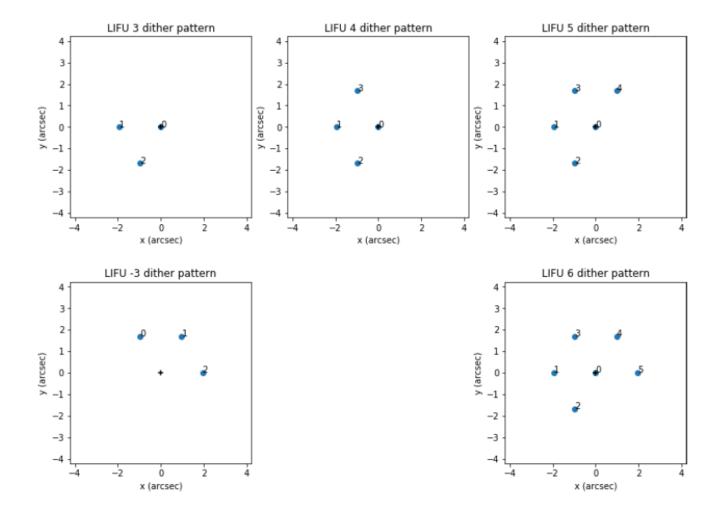

For custom dither patterns, constraints on the dither step size are imposed by the WASP to ensure that the guide star remains within the guiding camera field of view. From the perspective of IFU users, careful consideration of the dithering options should be made. Below it's an example of custom dithering. Note that pointings are provided for each dither, rather than the dither pointings from an initial position. The three pointings provided in the example below would be executed as part of one only OB.

```
# A custom 3-dither pattern LIFU pointing

TARGID,TARGNAME,GAIA_RA,GAIA_DEC,PROGTEMP,OBSTEMP,IFU_DITHER
PSZ1_125A,PSZ1_125A,316.369609537,-4.71060356792,41331,IAEEB,-1
PSZ1_125A,PSZ1_125A,316.369058752,-4.70967250147,41331,IAEEB,-1
PSZ1_125A,PSZ1_125A,316.370143237,-4.70966298315,41331,IAEEB,-1
```

Summary of table columns and examples


Column	Description	Format	Length	Value(s)	Units	Example
TARGPRO	Optional description of programme	ASCII	40			WS2023B2-010_001
TARGNAME	The target name	ASCII	30			M33
TARGID	The identifier of the target assigned for this programme	ASCII	30			M33_NE
PROGTEMP	Observing Programme Template	ASCII	8 (fixed)	http://wasp.ast.cam.ac.uk /progtemp		41331.3+
OBSTEMP	Observing Constraints Template	ASCII	5 (fixed)	http://wasp.ast.cam.ac.uk /obstemp		FAACA
GAIA_ID	GAIA Source Identifier	ASCII	19			14639761959685356 80


GAIA_DR	GAIA Data Release version	I1	2, 3		3
GAIA_RA	Gaia RA of target	F11.7	0360	degrees	178.221875
GAIA_DEC	Gaia Dec of target	F11.7	-9090	degrees	44.123919
GAIA_EPOCH	Gaia Epoch of target	F6.1	2015.5 (DR2), 2016.0 (DR3)	Julian year	2016.0
GAIA_PMRA	Gaia Proper Motion of target in RA	F11.3		mas/yr	12.1
GAIA_PMDEC	Gaia Proper Motion of target in Dec	F11.3		mas/yr	0.01
GAIA_PARAL	Gaia Parallax of target	F10.3		mas	0.002
IFU_PA_REQUEST	Position Angle of IFU bundle	F11.7	-180180	degrees	106.701
IFU_DITHER	IFU dither pattern code	12	-3, -1, 0, 3, 4, 5, 6		3


Aladin overlays of the LIFU footprints

Aladin overlays can be useful for preparing your observations. To use the Aladin overlays below, follow these instructions:

- Download the vot files linked below.
 Open Aladin and enter an object name in the command field, ex. 'M51'.
- 3. File Load instrument FoV.
- 4. In the emerging window, click on 'File' and browse for a suitable vot file.5. Click to submit. The vot file will load as a layer with a given colour.

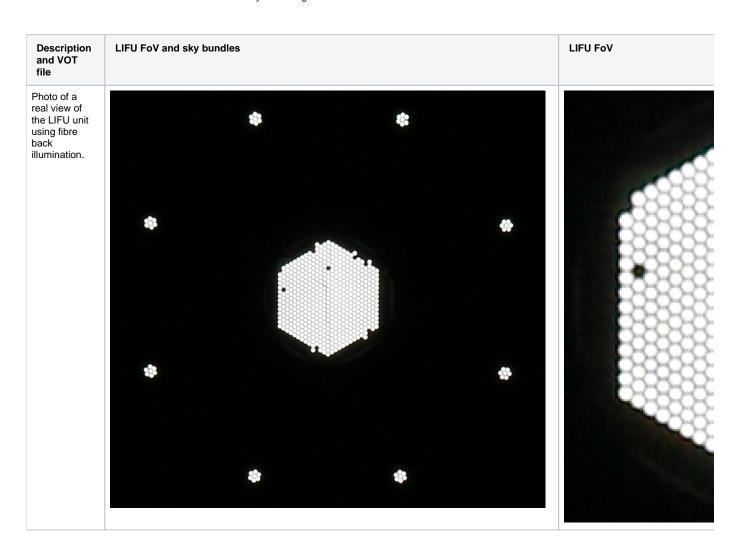

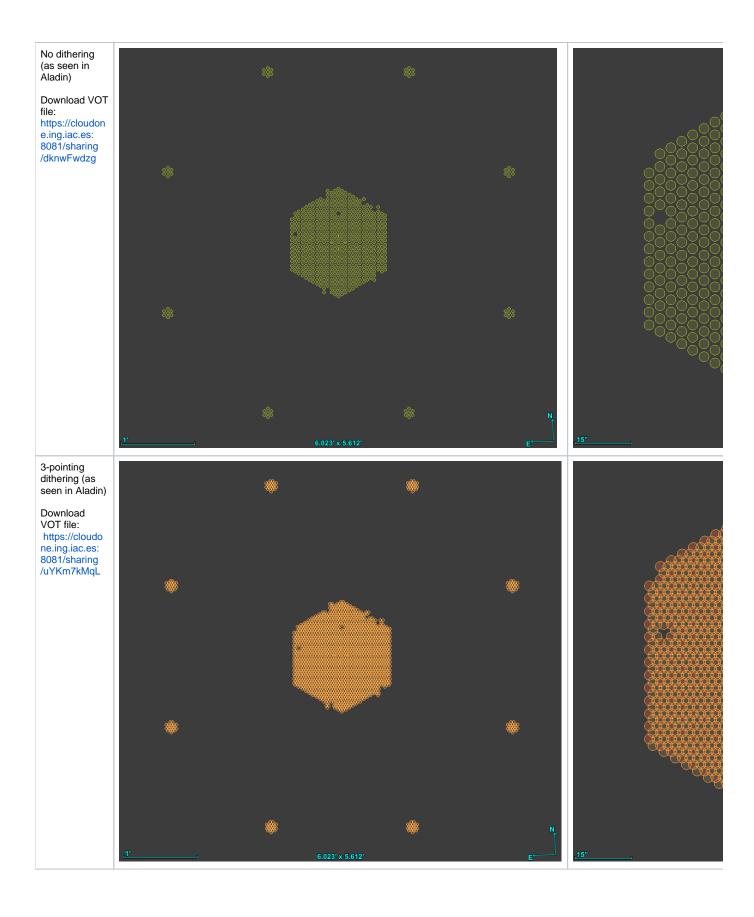

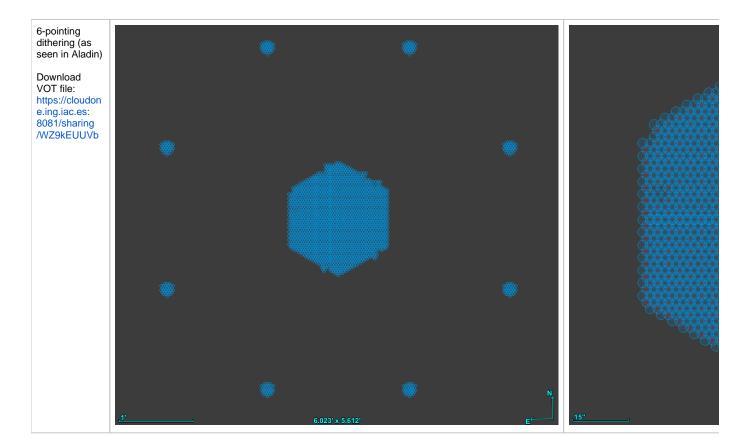

Catalogue upload to WASP

The edited catalogue must go through the IFU workflow, a set of scripts which will make a version of the FITS catalogue suitable for uploading to WASP (http://wasp.ast.cam.ac.uk/catsubmit/upload) for validation. First log on using the provided credentials (the ones emailed from WASP), note that you need to register your IP at http://wasp.ast.cam.ac.uk/ipsubmit/upload first. See below the tab used for uploading the catalogue.



At ING, we can provide assistance in guiding you through the different stages of the IFU workflow package, or we can even process your catalogue provided in any of the formats described above, on your behalf. In this latter case, simply email your catalogue to the contact below:

weave_open_time_support@ing.iac.es

LIFU mode: OB preparation

Download the validated catalogue

After the upload of the catalogue to WASP, this is validated. When it is passed, then the PI can download a FITS file (linked at the Download column in the figure below) which is the same catalogue as the one uploaded but this time including CNAMEs for every row entry.

Generation of OBs

Submission Management

The downloaded validated FITS catalogue is used by the IFU workflow software to generate OBs in the form of XML files. Every OB defines the observations for a given target, including instrument setup, observing conditions, dithering (if any), and exposure times. There are a number of considerations which don't affect the generation of an OB: if an OB requires being observed a number of times (as defined in the PROGTEMP code) then that OB will be cloned accordingly when ingested into the observing system. Also, if these OBs are chained: the observing system will handle the respective internal priorities.

At ING, we can provide assistance in generating the OBs from the CNAMEd catalogue, please contact:

Uploads to the WASP are managed in the Submission Management page. Please go to that page in order to view your past submissions and select those you wish to observe with

weave_open_time_support@ing.iac.es

OB upload to WASP

The processed OBs can be then uploaded, as a tar file if many, to http://wasp.ast.cam.ac.uk/xmlsubmit/upload (see below). Remember to log on first.



In the figure below we show an example of 6 OBs uploaded whose validation failed (the generated report will tell you why). The "Submission Management" button gives access you to all the catalogue and OB uploads made by the user. Note that if you have requested OB repetitions (for instance, PROGTEMP=41331.5, or 5 OB copies), then only one OB is validated on WASP (not 5 OBs following the example). The OB copies are actually made at the time they are ingested into the observing queue according to the information encoded in PROGTEMP.

IMPORTANT: Don't forget to select the right OBs in the submission maganement! Otherwise, they won't be queued.

JobID 1678287807

#	ID	File 6		User	Field	RA	Dec	PROGTEMP	OBSTEMP	Plate	File contents	Process State	Result	Report	Download
1	682	1678287807.tar	÷	ws2023a1-022		-					6 OBs	COMPLETED	FAILED	自±☎	
2	683	1678287807_1.xml	±	ws2023a1-022	L/022/Abell 2142/A2142center/d0	239.6067	27.1407	41332	NBDEB	LIFU	G1S56C0T547	COMPLETED	FAILED	ê± ≅	
3	684	1678287807_2.xml	÷	ws2023a1-022	L/022/Abell 2142/A2142center/d3	239.6067	27.1407	41332	NBDEB	LIFU	G1S56C0T547	COMPLETED	FAILED		
4	685	1678287807_3.xml	ŧ	ws2023a1-022	L/022/Abell 2142A/A2142infall/d0	239.7121	27.2988	41332	NBDEB	LIFU	G1S56C0T547	COMPLETED	FAILED		
5	686	1678287807_4.xml	ŧ	ws2023a1-022	L/022/Abell 2142A/A2142infall/d3	239.7121	27.2988	41332	NBDEB	LIFU	G1S56C0T547	COMPLETED	FAILED		
6	687	1678287807_5.xml	ŧ	ws2023a1-022	L/022/Abell 2142A/A2142infall/d6	239.7121	27.2988	61662	NBDEB	LIFU	G1S56C0T547	COMPLETED	FAILED	自±⊠	
7	688	1678287807_6.xml	ŧ	ws2023a1-022	m/022/CLUSTER1/d3	316.2583	-5.0368	71331	IAEEB	mIFU	G7S104C1T37	COMPLETED	FAILED	ê ± ≅	

Uploads to the WASP are managed in the Submission Management page. Please go to that page in order to view your past submissions and select those you wish to observe with.

Submission Management

Observations and follow up

ING ingests the OBs (as provided by CASU on WASP) into the observing system at the start of every observing period. OBs are then queued taking into account the priorities derived from the band the associated programme lies (A, B or C). See the open-time time allocations. Observations can be followed up at CASU's operational repository.